

- Two configurations available to fit each patient's incontinence condition
- Minimally invasive adjustment system



The constitution of the urethra and the pelvic floor response varies from patient to patient.

Artificial Urinary Sphincters (AUS) have been established as the Gold Standard in the surgical treatment of male stress urinary incontinence. While former concepts of the AUS have failed to adapt to the variety of urethras and pelvic floor conditions, we understand that each patient deserves the most individual treatment possible. As a result, VICTO Adjustable Artificial Sphincter is available in two configurations:





VICTO and VICTO+ are both adjustable artificial urinary sphincters designed to treat male stress urinary incontinence.

VICTO technology and features make it the most evolved AUS on the market:

MINIMALLY INVASIVE ADJUSTMENT SYSTEM

Both VICTO configurations come with an adjustment system that allows for unlimited easy-to-perform pressure adjustments of the already implanted product. This procedure can be performed on outpatient basis.

UNIQUE PUMP DESIGN

The VICTO pump, made of soft silicone, has an easy-to-palpate design which makes it very easy for the patient to operate.

PRE-CONNECTED DEVICE

VICTO and VICTO+ are provided pre-connected, ready for filling and implantation.

CIRCUMFERENTIAL DESIGN

Homogeneous pressure around the urethra

MULTIPLE CUFF SIZES

To maximize product performance while minimizing the risk of erosion, a variery of cuff size options that provide uniform pressure around the urethra are required.

VICTO offers four cuff sizes which can be chosen according to the anatomy of each patient.


















High and permanent occlusive pressure, especially on severely damaged pelvic floor structures, is considered to be the main reason for complications over time.

Patients with a dysfunctional pelvic floor do not have the ability to disrupt their urinary flow. These patients need additional pressure in the cuff to keep the urethra closed in physical stress situations.

VICTO+ is the perfect solution for these patients because it comes with an extra balloon that provides increased cuff pressure on demand by forwarding abdominal pressure to the cuff in stress conditions.

Working together, the easy adjustment system and the unique dynamic response feature in VICTO+ allow a lower baseline pressure while at the same time keeping patients with a dysfunctional pelvic floor continent in stress situations.

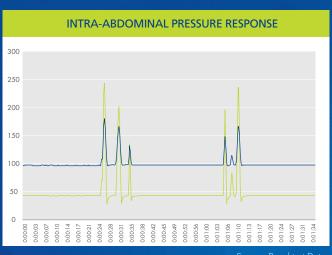

Source: Benchtest Data



ORDERING INFORMATION

Starter Kits*		
DESCRIPTION	CUFF SIZE	CODE
VICTO STARTER KIT 3.7	3.7	KIT-AUS-37-WS
VICTO STARTER KIT 4.0	4.0	KIT-AUS-40-WS
VICTO STARTER KIT 4.5	4.5	KIT-AUS-45-WS
VICTO STARTER KIT 5.0	5.0	KIT-AUS-50-WS
VICTO+ STARTER KIT 3.7	3.7	KIT-AUS-37
VICTO+ STARTER KIT 4.0	4.0	KIT-AUS-40
VICTO+ STARTER KIT 4.5	4.5	KIT-AUS-45
VICTO+ STARTER KIT 5.0	5.0	KIT-AUS-50

^{*}Each Starter Kit contains an empty AUS system in the respective size and comes with the needles and syringes to fill the system.

Accessories Kit	
DESCRIPTION	CODE
VICTO / VICTO+ ACCESSORIES KIT	KIT-AC-01

Replacement Kit	
PUMPS	CODE
VICTO / VICTO+ PUMP	KIT-RP-03-10-06
BALLOONS	CODE
VICTO / VICTO+ BALLOON	KIT-RB-O1
CUFFS	CODE
VICTO / VICTO+ CUFF 3.7	KIT-RC-37-10-06
VICTO / VICTO+ CUFF 4.0	KIT-RC-40-10-06
VICTO / VICTO+ CUFF 4.5	KIT-RC-45-10-06
VICTO / VICTO+ CUFF 5.0	KIT-RC-50-10-06

Adjustment Kit	
DESCRIPTION	CODE
VICTO / VICTO+ ADJUSTMENT KIT	KIT-AD-06

victosphincter.com

